

HandyStep[®] touch

The next generation repetitive pipette

Simple

Touchscreen operation

Easy

Automatic tip ejection

Efficient

Compatible with various third party tips



BRAND. For lab. For life.

BRANDTECH[®]
SCIENTIFIC, INC.

HandyStep® touch · HandyStep® touch S

Overview

- Touchscreen operation with intuitive menu navigation and integrated help function
- Two models available: HandyStep® touch or HandyStep® touch S with expanded functionality
- Menu languages: English, German, Spanish, French and Chinese
- Variable volume adjustment from 1.0 µL to 50 mL
- Inductive charging with optional charging stand

A closer look

Versatile

HandyStep® touch with the standard modes Multi-Dispensing, Auto-Dispensing, Pipetting and Favorites;

HandyStep® touch S with the additional functions Sequential Dispensing, Multi-Aspiration and Titration

Long-lasting

Many hours of continuous operation with a single battery charge

Robust

High-quality materials provide long service life, and smooth surfaces make cleaning/ disinfection easy



USB-C connection

For quick charging by cable

STEP button

Ergonomically arranged STEP button for effortless volume filling and delivery



Easy

Automatic size detection of BRAND® PD-Tip™ II

Efficient

In addition to BRAND® PD-Tip™ II, compatible third party tips can also be used.

Simple

Automatic tip ejection by the press of a button prevents contamination



Convenient

Inductive charging with the optional charging stand

Charging stand with HandyStep® touch S and PD-Tip™ II

Intuitive operation

Be flexible and have the right solution at hand for every application. Even complex procedures are easy to carry out – no problem with the expanded functionality of the HandyStep® touch S.

Mode	Function	HandyStep® touch	HandyStep® touch S
Multi Dispensing (MULTI-DISP)	Successive dispensing of a set volume	✓	✓
Auto Dispensing (AUTO-DISP)	Automatic dispensing by the instrument	✓	✓
Pipetting (PIP)	Pipetting a defined volume	✓	✓
Sequential Dispensing (SEQ-DISP)	Successive dispensing of different volumes	–	✓
Multi Aspirating (MULTI-ASP)	Aspirating different volumes	–	✓
Titrating	Manual or successive delivery of liquids for titration	–	✓
Favorites	Saving favorites for quick access	✓	✓
Settings	Individual customization of the device (e.g. brightness or device name)	✓	✓



Expanded functionality on the HandyStep® touch S:

The clearly structured menus make fast work with intuitive operation possible, while the most important functions and parameters are always in view.




Easy to understand

Main menu with all modes for quick selection



Clearly structured

Working mode with all important information

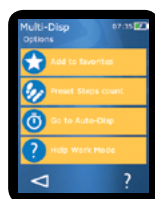

Organized

All options visible at a glance



Helpful

Internal help right where it is needed



Multi dispensing with the HandyStep® touch and PD-Tip™ //

The perfect team

Automatic size detection with PD-Tip™ // from BRAND®

The HandyStep® touch saves time and prevents errors through automatic tip size recognition of the PD-Tip™ // from BRAND®, with patented coding on the pistons. After inserting the tip, the size is automatically recognized and displayed, making it easy to select the volume to be dispensed. When a new PD-Tip™ of the same size is inserted, all instrument settings are maintained.



Sequential dispensing with the HandyStep® touch S and PD-Tip™ //

Technical data

HandyStep® touch HandyStep® touch S

Repetitive pipette, includes universal AC adapter (for USA, Europe, UK, Japan, Australia and New Zealand), shelf/rack mount, sample PD-Tip™ II, operating manual and performance certificate.



Support stand with HandyStep® touch



HandyStep® touch

HandyStep® touch S

Description	Cat. No.	2019 List Price
HandyStep® touch	705200	\$1,055.00
HandyStep® touch S	705210	1,299.00

Accessories

Description	Cat. No.	2019 List Price
Charging stand	705220	\$345.00
Support stand	705230	144.00
Shelf/rack mount	705235	35.00
Universal AC adapter	705250	85.00
Spare battery	705225	80.00

PD-Tip™ II

Precision dispenser tips coded for automatic size recognition.

Capacity mL	Pack of	non-sterile Cat. No.	2019 List Price	BIO-CERT® Cat. No.	2019 List Price
0.1	100	705700	\$124.00	705730	\$170.00
0.5	100	705702	124.00	705732	170.00
1.0	100	705704	124.00	705734	170.00
1.25	100	705706	124.00	705736	170.00
2.5	100	705708	124.00	705738	170.00
5	100	705710	124.00	705740	170.00
10	100	705712	124.00	705742	170.00
12.5	100	705714	124.00	705744	170.00
25	50 + 1 adapter	705716	84.00	705746*	52.00
50	25 + 1 adapter	705718	42.00	705748	52.00
Set of PD-Tip™ II, non-sterile, 20 pieces per size: 0.5/1/1.25/2.5/5/10/12.5 mL				705720	168.00
Adapter for 25 and 50 mL PD-Tip™ & PD-Tip™ II, PP, autoclavable		Pack of 10 702398	96.80	Pack of 5 702399	76.20

BIO-CERT® : Pack of 25 + 1 adapter



For more information visit www.brandtech.com or scan the QR Code.

HandyStep® touch, HandyStep® touch S and PD-Tips II have patents & patent applied for features. California Residents: For more information concerning California Proposition 65, please refer to www.brand.de/calprop65.

Subject to technical modification without notice.



Lab Rats Trust BrandTech!

BRANDTECH®
SCIENTIFIC, INC.

888-522-2726
www.brandtech.com